


PRIMACORE™

BIOPSY GUN - SEMI AUTOMATIC

Designed for **smooth** and **efficient** collection of **exact, high quality** core samples from soft tissues and organs.

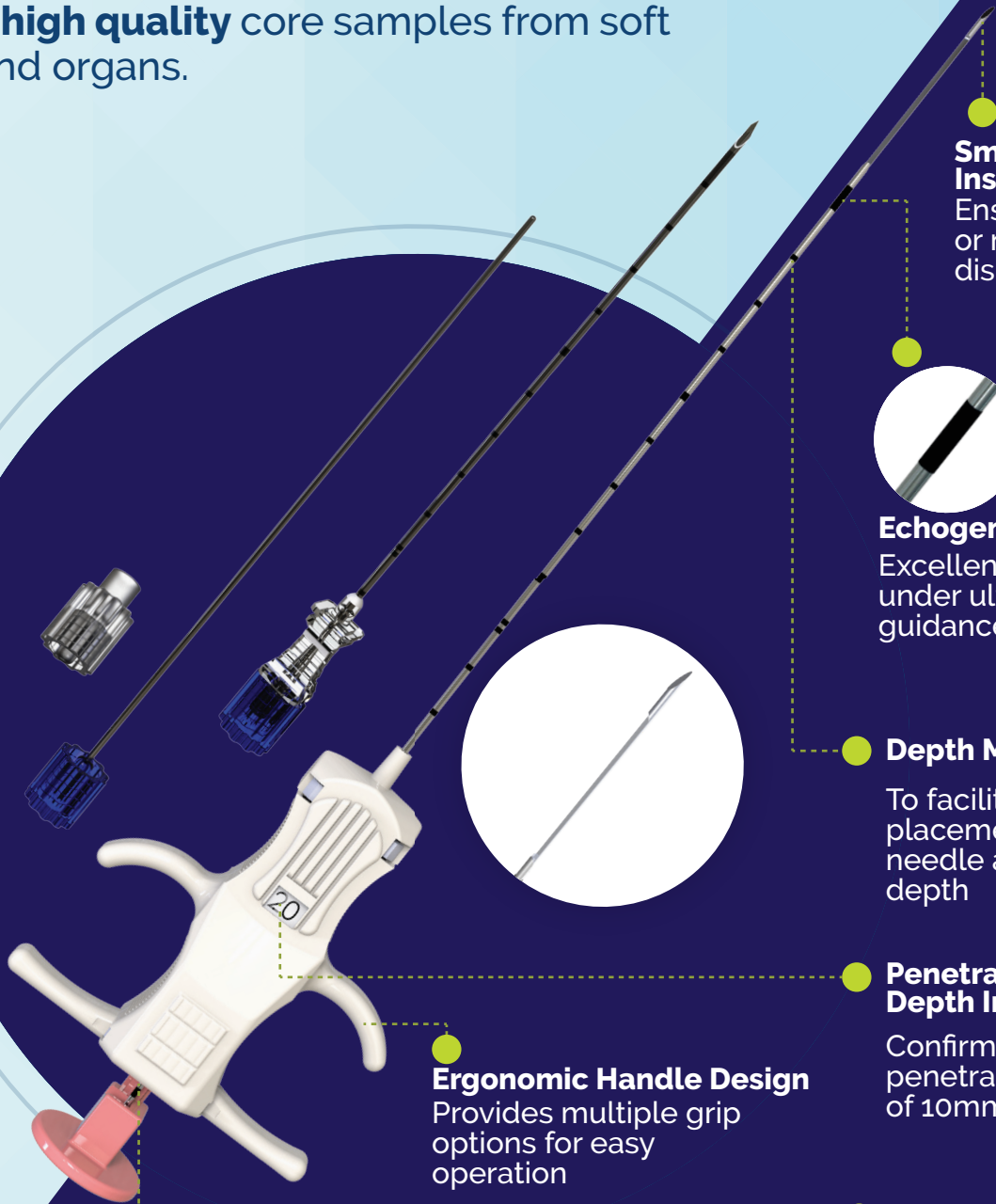

Smooth Stylet Insertion
Ensures minimal or no tissue displacement



Echogenic Needle Tip
Excellent tip visibility under ultrasound guidance

Depth Markings
To facilitate the precise placement of the needle at the desired depth

Penetration Depth Indicator
Confirms primed penetration depth of 10mm or 20mm

Ergonomic Handle Design
Provides multiple grip options for easy operation

Fire-ready Indicator
Designed to reduce the risk of premature firing

SMOOTH access:

- Ergonomic handle with multiple grip options
- Excellent tip geometry to reduce penetration force and tissue displacement
- Exemplary tip visibility under ultrasound guidance

EFFICIENT operation:

- Convenient adjustable throw lengths – 10mm and 20mm
- Smooth stylet insertion
- Siliconized needles to reduce needle deflection
- High Stylet stiffness avoids needle bending

EXACT, High-quality sample collection:

- Penetration depth indicator – confirms primed penetration depth
- Semi-automatic firing – aids visual verification of sample notch placement
- Superlative tactile feedback of tissue penetration and collection

Complete Instrument Kit:

- Primacore™ Biopsy Instrument
- Compatible Coaxial
- Trocar and Blunt Tip Stylet
- Depth Stopper
- 10 mm Adapter

Product Specifications

| Size (G) | Length (cm) |
|----------|-----------------|
| 14 | 10, 16 & 20 |
| 16 | 10, 16 & 20 |
| 18 | 10, 16, 20 & 25 |
| 20 | 10, 16 & 20 |

BLUENEEM®
CLINICIANS' TRUSTED PARTNER



Blue Neem Medical Devices Pvt Ltd
Harohalli Industrial Area, Bengaluru - 562 112, Karnataka, India



+91 97399 72855



enquiry@blueneem.com



www.blueneem.com



360° Product
Experience