

DALELA TRIPLE J URETERAL STENT

(PATENT PENDING**)

Key Features

- An additional 'J' like bend is created at bladder end, as a result of this bend the lower coil of stent gets directed towards the dome and away from the sensitive trigone.
- Thermosensitive, advanced PU material used for construction
 - Stent will be stiffer at room temperature promoting pushability for stent deployment
 - Material becomes softer at body temperature enhancing patient comfort

Patient Benefits


- The design modification and the material features, reduces dysuria & other stent induced symptoms
- Enhances patient comfort
- Long-term Indwelling
- Radiopaque for better visualization
- Easy removal as the lower coil seen easily towards dome

Size (in Fr)

3, 3.5, 3.7, 3.8, 4, 4.5, 4.7, 4.8, 5, 5.5, 6, 6.5
7 & 8

Length (in cm)

22, 24 & 26

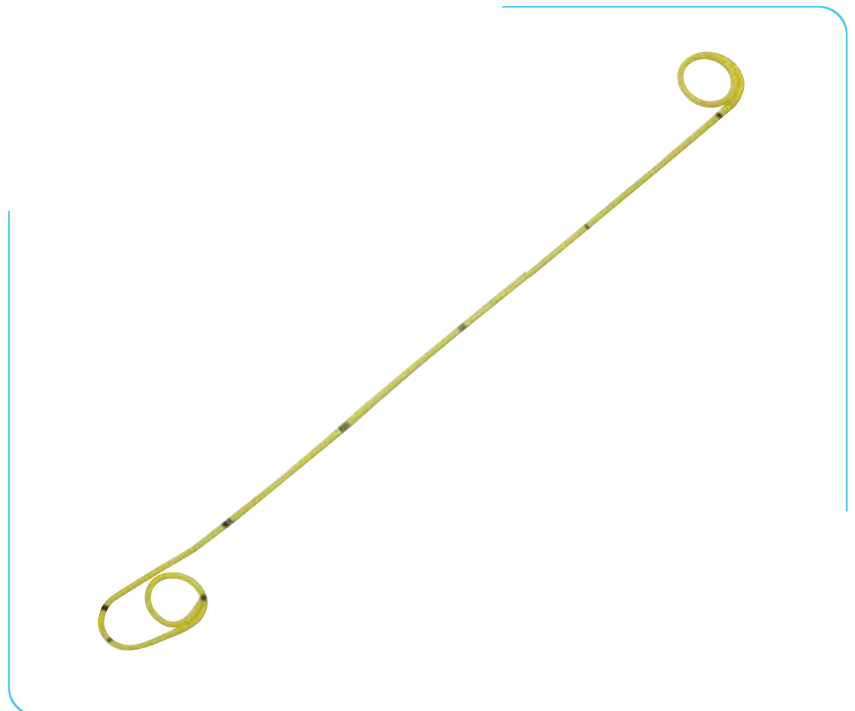

Note: During stent placement, deploy the lower end in full bladder to let the lower end go towards the dome.