


## Kidsoft<sup>®</sup> Sigma Catheter Kit

### Intended Use:

Used to drain urine, internal or external, after pyeloplasty procedure.

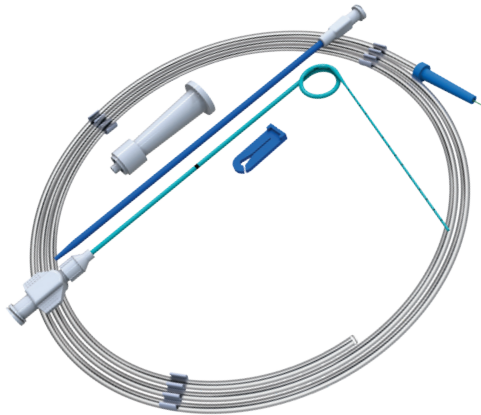

| Size (in Fr)    | Length (in cm) |
|-----------------|----------------|
| 2.8, 3, 3.5 & 4 | 10/15          |

### Set Includes:

PTFE Guidewire (Straight), Catheter, Dilator, Connector Hub & Pinch Clamp

### Key Features and Benefits:

- The central 'pigtail' remains in the pelvis while the narrow distal extension stents the anastomosis and the ureter
- The Distal Ureteric end can be altered/trimmed as per the requirement.
- It also has centimetre Marking to help healthcare providers to ensure correct placement or insertion depth.
- 90 degree angle avoids stent migration
- The proximal extension tubing can be used for external drainage of urine from the kidney or can be blocked by using the clamp provided
- Prevents post-operative bladder spasms
- Can be removed on an out-patient basis
- Avoids subsequent surgical intervention for catheter removal with its attendant costs and potential for infections
- Extension dia 6Fr