


KIDSOFT LAPAROSCOPIC PYELOPLASTY STENT (URETERAL STENT)

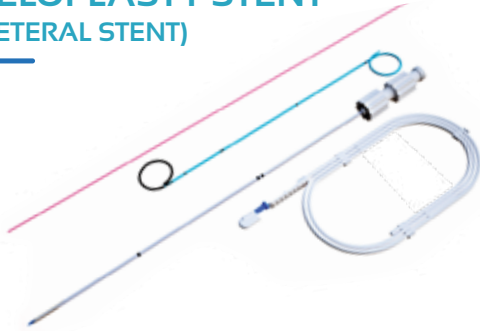

Set Includes

Introducer with Needle,
Guidewire, Stent
& Pusher

Size (in Fr)

2.8, 3, 3.5, 4, 5 & 6

Length (in cm)

12, 14, 15 & 16

Intended Use:

Intended to drain urine after pyeloplasty procedure.

Key Features:

- Ergonomically designed for stent placement during Laparoscopic Pyeloplasty
- The needle can be percutaneously introduced and placed across the anastomosis after removing the sharp stylet
- The stent can be guided into the ureter with ease, the Proximal end can be brought out as a nephrostomy if introduced through the Kidney or Pelvis
- The stent can be retained for upto 30 days