


R"EA"Plogel Double Lumen Esophagus Suction Tube

Intended Use:

A double lumen “sump” suction tube is used in the management of neonates and infants with oesophageal atresia and is used to continuously suck the saliva from the upper oesophageal pouch without causing mucosal damage.

Will be very useful during transportation and for pre-operative management of babies with Esophageal Atresia.

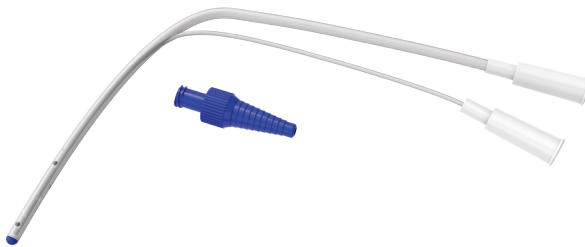

Package Includes
R “EA” Plogel Tube and Catheter Adaptor

Size (in Fr)
12/6

Length (in cm)
25

Key Features and Benefits:

- Clear dual-lumen tube with suction adaptor for both the lumens
- Rounded atraumatic tip
- Wide side holes for ease of suction without mucosal damage
- Multi holes at the distal tip for enhanced drainage