



WORTX
URETERAL ACCESS SHEATH
WITH SUCTION

Suction.
Evacuation.
Collection.

WortX™ facilitates ureteral dilation, ensuring smooth and efficient scope access to the ureter. It enhances safety and effectiveness during RIRS procedures, contributing to an improved stone-free rate (SFR).

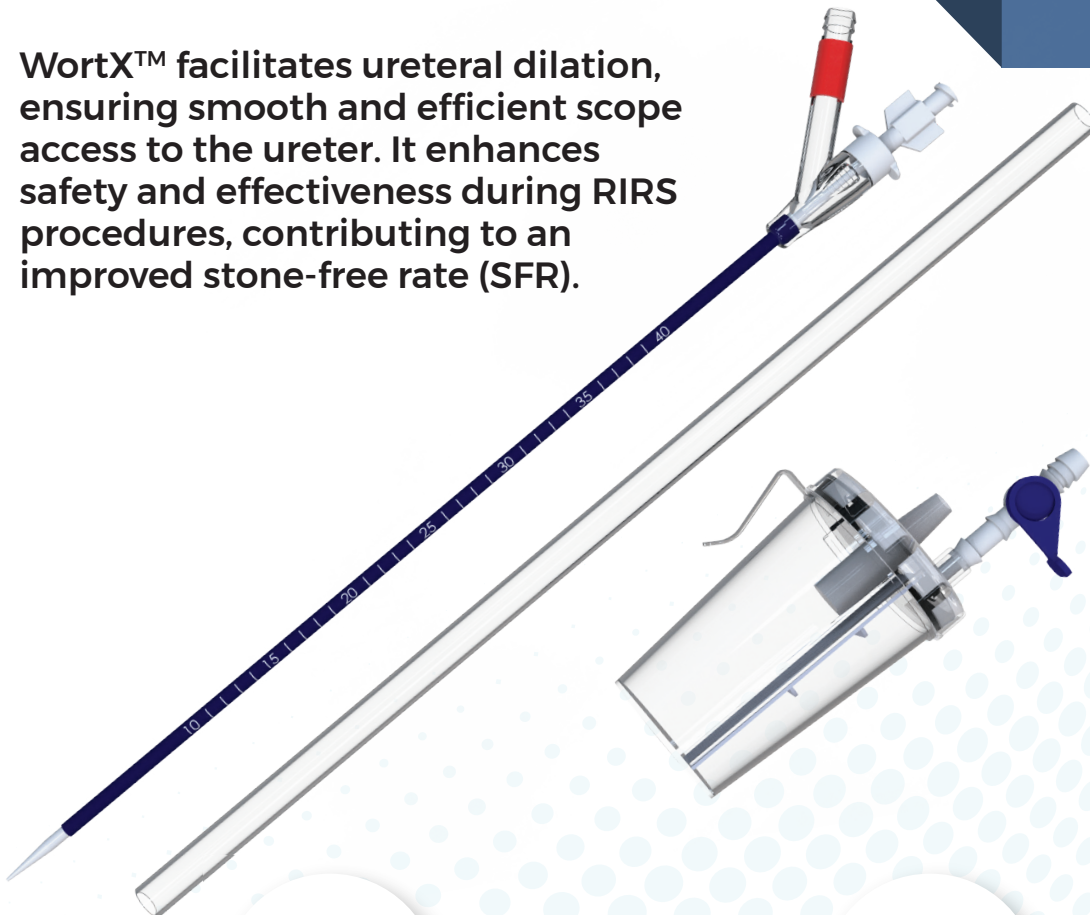

Flexible sheath
for easy access

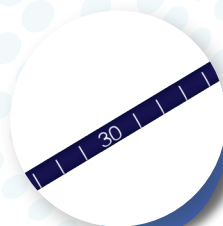

Length markings
for accurate
placement



Hydrophilic coating
for easy insertion



Tapered tip for
smoother insertion




Key Features And Benefits



Flexible Sheath

10cm flexible tip for easy kidney access with a scope, a rigid shaft for scope support

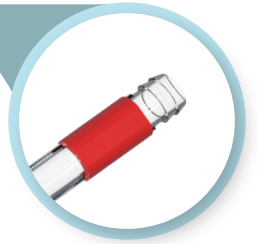

Suction Vent

Integrated suction continuously removes fluids and debris, enhancing visibility and shortening procedure time



Air-Tight Seal

The bottle features an airtight seal to ensure effective suction throughout the procedure



Graduated Markings

Markings on the side measure the volume of the fluid



Collector

Collects debris and prevents it from entering the suction machine

Specifications

1. OD of the UAS - 12 Fr
2. ID of the UAS - 10 Fr
3. Length of the UAS - 50 cm
4. Use Guidewire Size - 0.035 Inch
5. Stone Collection Bottle - 200 ml
6. Silicon Tube - Secure Fit for Stone Collection Bottle and Device

BLUENEEM[®]
CLINICIANS' TRUSTED PARTNER



Blue Neem Medical Devices Pvt Ltd
Harohalli Industrial Area
Bengaluru - 562 112, Karnataka, India



+91 97399 72855



enquiry@blueneem.com