



**Suction.  
Evacuation.  
Collection.**

WortX™ facilitates ureteral dilation, ensuring smooth and efficient scope access to the ureter. It enhances safety and effectiveness during RIRS procedures, contributing to an improved stone-free rate (SFR).

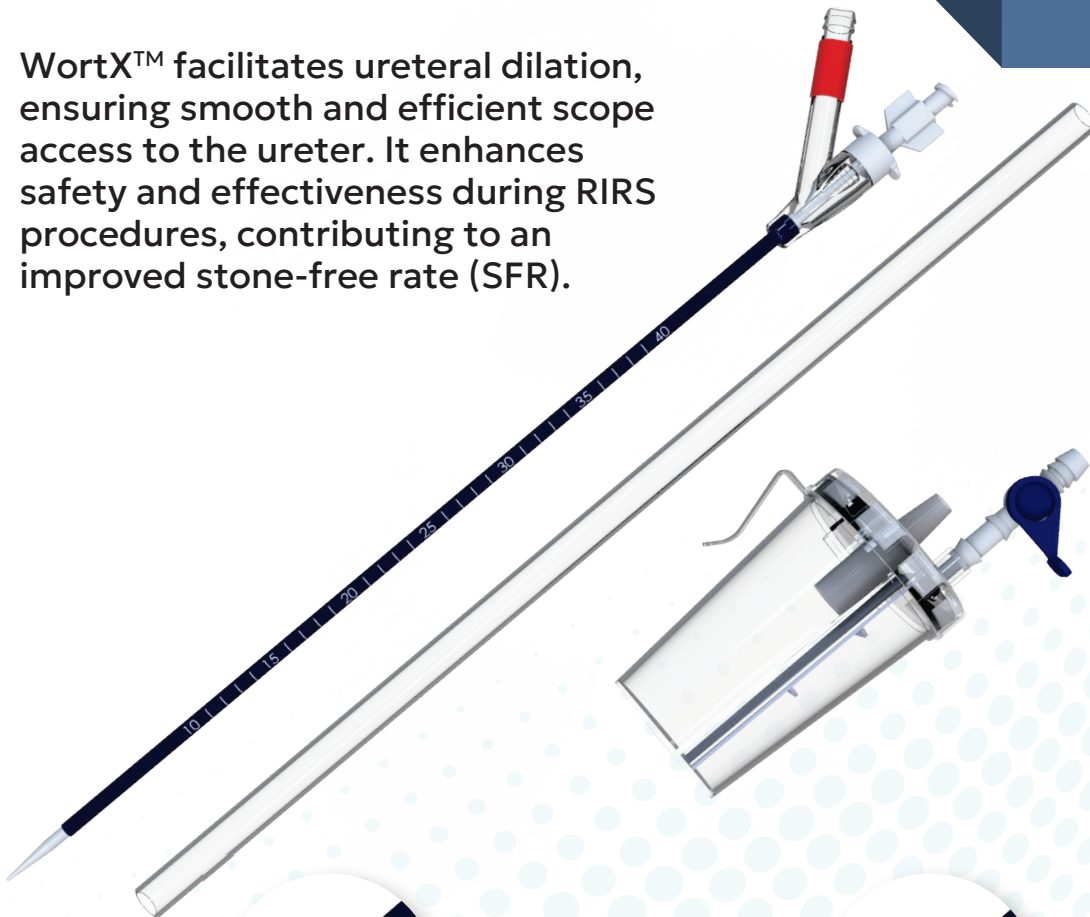

**Flexible sheath  
for easy access**

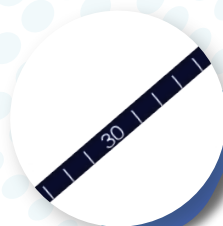

**Length markings  
for accurate  
placement**



**Hydrophilic coating  
for easy insertion**



**Tapered tip for  
smoother insertion**




## Key Features And Benefits



### Flexible Sheath

10cm flexible tip for easy kidney access with a scope, a rigid shaft for scope support

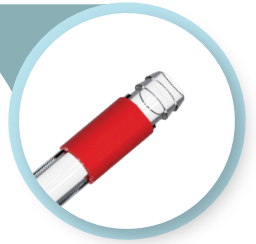


### Suction Vent

Integrated suction continuously removes fluids and debris, enhancing visibility and shortening procedure time



### Air-Tight Seal

The bottle features an airtight seal to ensure effective suction throughout the procedure

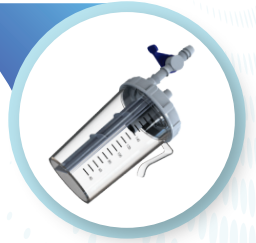

### Graduated Markings

Markings on the side measure the volume of the fluid



### Collector

Collects debris and prevents it from entering the suction machine

## Set Includes

1. OD of the UAS – 12 Fr
2. ID of the UAS – 10 Fr
3. Length of the UAS – 35, 40, 45 & 50 cm
4. Use Guidewire Size – 0.035 Inch
5. Stone Collection Bottle – 200 ml
6. Silicon Tube – Secure Fit for Stone Collection Bottle and Device

**BLUENEEM**<sup>®</sup>  
CLINICIANS' TRUSTED PARTNER



**Blue Neem Medical Devices Pvt Ltd**  
Plot Nos 270 & 271, Road No 5  
Harohalli Industrial Area, Phase II  
Kanakapura Taluk  
Ramanagara Dist - 562112  
Karnataka, India



+91 97399 72855



[enquiry@blueneem.com](mailto:enquiry@blueneem.com)