




Intended Use:

Intended to be implanted for drainage of urine from renal pelvis to urinary bladder

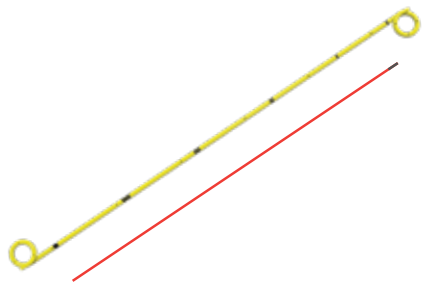

Types

Both Ends Open, One End Closed, Both Ends Multiloop, One End Multiloop
With or Without Hydrophilic coating, With or Without String

Package Includes
Stent with pusher

Also available in set, which Includes
Stent with pusher and guidewire

Key Features and Benefits:

- Thermosensitive, advanced PU material for enhanced patient comfort
 - Stent will be stiffer at room temperature promoting pushability for stent deployment
 - Material becomes softer at body temperature
- Optimal loop shape and high memory retention reduces the risk of migration
- Optimal radiopacity for better visualization

Type	Size (in Fr)	Length (in cm)
Standard	4.5, 4.7, 4.8, 5, 5.5, 6 6.5, 7 & 8	12, 14, 16, 18, 20, 22, 24 26, 28 & 30
Multiloop	4.7, 4.8, 5, 6, 7 & 8	18-24 & 24-30
Pediatric	2.8, 3, 3.5, 3.7, 3.8 & 4	6, 8, 10, 12, 14 & 16