


PCN Catheter With Needle

Intended Use:

Used for percutaneous placement of a pigtail catheter in the renal pelvis for nephrostomy drainage by direct stick technique.

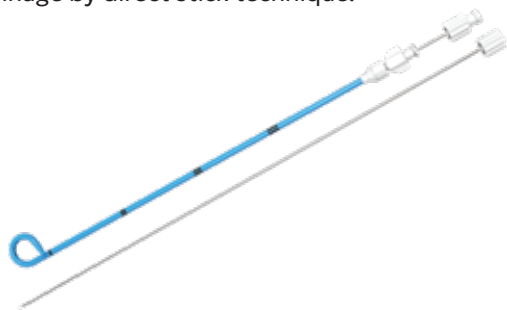

Package Includes

PCN Catheter with metal straightener, metal trocar and urine bag connector

Key Features and Benefits:

- The pigtail tip configuration helps in the retention of the distal tip in the renal pelvis
- Eye holes for better drainage
- Trocar design and catheter tip design prevents buckling of catheter during the procedure
- Also available in Hydrophilic coating for easy insertion
- Also available with locking mechanism

Type	Size (in Fr)	Length (in cm)
PCN Catheter with Needle	6, 7, 8, 8.5, 9, 10, 12, 14 16 & 18	20, 22, 25, 30 & 35
PCN Catheter with Needle and Locking Mechanism	6, 7, 8, 10, 12, 14 & 16	20, 22, 25, 30 & 35